
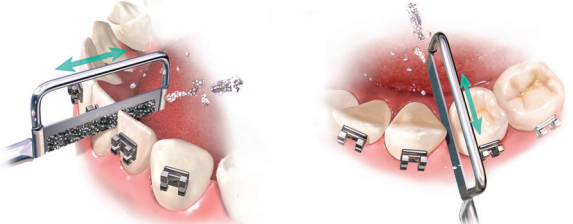






SPECIALTY HEADS – For Specialty Needs

Interproximal Stripping and Finishing

Clinical Application	Special Movement	Head	
For interproximal polishing, stripping or preparation with Ortho-Strips or Proxoshape Diamond Files (from Axis Dental).	With attachments 20CH and 20LP 20,000 rpm	 <p>61LRG 36 individual locking positions File elevation of 0.8mm Buffered stroke (Protects soft tissue) Order #0549.3450</p>	 
For interproximal cleaning, polishing, and preparation with Proxoshape Diamond Files (from Axis Dental).	With attachments 20CH and 20LP 20,000 rpm	 <p>61LG Free rotation of files File elevation of 0.8mm Buffered stroke (Protects soft tissue) Order #0549.3410</p>	 

A
"MUST HAVE"
if you offer
INVISALIGN®*
to your patients

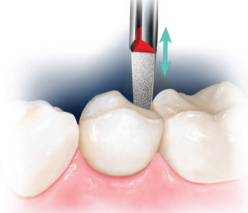

Ortho-Strips

Intensiv Ortho-Strips from AXIS Dental are the solution for one of the most challenging aspects of orthodontic procedures, interproximal stripping. Intensiv Ortho-Strips, together with the 61LRG low-speed stripping. Intensiv Ortho-Strips, together with the 61LRG low-speed head from KaVo, provide a predictable solution to this problematic procedure. Available in five grits for initial interproximal stripping to interproximal finishing.

Proxoshape Diamond Files

Intensive Proxoshape Diamond Files from AXIS Dental** can be used to contour and polish the enamel as desired after using the Intensiv Ortho-Strip system for interproximal stripping. Intensiv Proxoshape Diamond Files are also great for the removal of amalgam or composite overhangs. Available in numerous grits and lengths.

Specialty Endo Heads

Clinical Application	Speed at Bur up to	Head	Attachment/Motor
Mechanical root canal treatment with hand instruments (shaft diameter 3.5 mm to 4.0 mm)	Oscillating 90° movement, operates with attachments 7 and 29	 <p>3LD Push Button Endo head for hand files Order #0549.3720</p>	These heads fit any KaVo low-speed system.
Mechanical root canal treatment with latch instruments (shaft diameter 2.35 mm)	Oscillating 90° movement, operates with attachments 7 and 29	 <p>53LDN Push Button Endo head for latch type files Order #0549.3780</p>	

* Invisalign is a registered trademark of Align Technology.

** To reach an AXIS Dental product specialist call toll free 888-654-2947.